




Description:

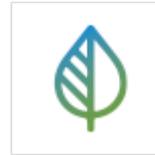

The Plus Series Pharmacy Freestanding Undercounter Freezers feature a microprocessor temperature controller and display, allow for precise temperature control, verification and recovery versus mechanical style thermostats. Critical samples and supplies are security backed with an extensive package of audible and visual alarms **Periodic defrosting is required.** **The Freestanding models ("UCFS") require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure.** **Environmentally friendly, natural hydrocarbon (HC) refrigerants (not containing hydrofluorocarbons) vastly reduce global-warming potential (GWP) while meeting new EPA/SNAP mandates and UL, ASHRAE compliance guidelines.** For more information on Vaccine Storage: [Vaccine Storage](#) *The acceptable temperature ranges for vaccine storage are 36°F to 46°F (or 2°C to 8°C) in the refrigerator and -58°F to 5°F (or -50°C and -15°C) in the freezer. **Adjustable operating temperature range: -4°F to -22°F (-20°C to -30°C). Factory Preset at -30°C. See "View/Print Product Information" drop down for instructions.***

Details

- THREE ADJUSTABLE SHELVES WITH GUARD RAIL ON BACK TO PREVENT PRODUCT FROM BEING PUSHED AGAINST THE BACK WALL
- BOTTOM FOOTED SHELF TO ALLOW AIR CIRCULATION
- BACK WALL EVAPORATOR COVER
- KEYED DOOR LOCK
- MAGNETIC DOOR GASKET FOR POSITIVE SEAL
- RECESSED DOOR HANDLE
- LEVELING LEGS
- LED INTERIOR LIGHT
- PHARMACY REFRIGERATOR/FREEZER TOOLKIT AND TEMPERATURE LOGS
- HFC-FREE REFRIGERANT (R290), FOAM INSULATION AND PACKAGING
- TWO YEAR PARTS AND LABOR WARRANTY, PLUS AN ADDITIONAL THREE YEAR COMPRESSOR PARTS ONLY WARRANTY
- **NOTE: THIS UNIT MUST HAVE 4" CLEARANCE ON SIDES AND BACK FOR ADEQUATE VENTILATION**

Specifications

- EXTERIOR DIMENSIONS: WIDTH: 23 ¾" DEPTH: 24" HEIGHT: 32 1/8"
- MANUAL DEFROST
- 115V, 60 HZ, 2.2 AMPS, 1/4 HP
- UL/C-UL LISTED
- SHIPPING WEIGHT: 128 LBS.
- **NOTE: THIS UNIT MUST HAVE 4" CLEARANCE ON SIDES AND BACK FOR ADEQUATE VENTILATION**

Temp Controller

- MICROPROCESSOR TEMPERATURE CONTROLLER
- DIGITAL TEMPERATURE DISPLAY AND ALARM MODULE WITH BATTERY BACK-UP TO SUPPORT TEMPERATURE MONITORING AND ALARMS
- AUDIBLE AND VISUAL HIGH AND LOW TEMPERATURE ALARMS
- REMOTE ALARM CONTACTS
- 2 PROBES; 1 IN AIR AND 1 IN GLYCOL BOTTLE
- ADJUSTABLE OPERATING TEMPERATURE RANGE: -4°F to -22°F (-20°C to -30°C) Factory Preset at -30°C
- Calibrated using a NIST traceable device, certificate included